

Keysight 5500 or 5100

MAC Mode[®] AFM Studies of Zinc-Induced DNA Kinking

Application Note

Introduction

The bending of DNA strands is a common phenomenon either intrinsic or induced by protein binding and/or other changes in the physiological environment (such as local ion concentrations). It is believed to play an important role in the control of gene expression, replication, recombination, and packaging in the nucleus. Kinking, roughly defined as abrupt bending, is a distinct category of DNA bending and is often observed in protein-DNA complexes. Numerous methods are available to detect structural fluctuations in DNA. To clearly distinguish between smoothly distributed bending and highly localized kinking, however, a spatial resolution on the order of at least 1 nm is required.

Among many high-resolution imaging tools, the atomic force microscope (AFM), with easy sample preparation and close to native imaging conditions, stands out as one of the best options for biologists. Its near-atomic resolution makes it useful for *in situ* imaging of conformational changes in DNA molecules. Due to tip-broadening effects, however, the width of DNA molecules obtained with atomic force microscopy has been on the order of 5 nm.

In this study, we employed microscopy, patented by Keysight Technologies, Inc., MAC Mode (Figure 1) atomic force, which uses a magnetically-driven oscillating probe with an oscillation amplitude significantly smaller than that of fluid tapping mode, or Tapping Mode. MAC Mode (Figures 2 and 3) was shown to image DNA molecules at a much higher resolution (about 1 nm) than other AFM techniques. We utilized MAC Mode to investigate DNA minicircles constructed from phased sequence motifs known to bend DNA.

Linear and circular molecules were separated with two-dimensional gels. A 168-bp DNA circle containing four of the above fragments was used for AFM study. DNA samples were diluted to a final concentration of 0.5 µg/ml using deionized water. Before imaging, MgCl₂ and/or ZnBr₂ were added to the DNA solution and the mixture was deposited on freshly cleaved mica of a Keysight AFM.

In a previous study using the same sample and technique, the 168-bp circles were shown to obtain dramatically different conformations in different cations. When imaged in 1 mM MgCl₂, the DNA circles were smoothly bent, all nearly circular. In 1 mM ZnBr₂, about 4 kinks on the average were observed in each molecule. In this work, we were mainly interested in the cation concentration-dependence of the DNA kinking phenomenon. In particular, what is the Zn²⁺ concentration threshold that induces kinking?

To answer this question, ZnBr₂ concentration was gradually decreased from 1 to 0.05 mM in a series of experiments.

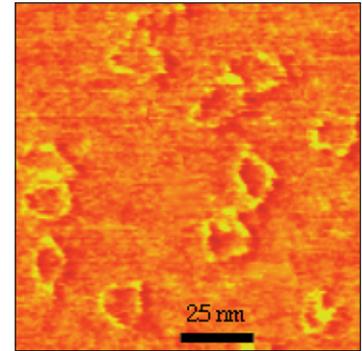


Figure 1. MAC Mode AFM image of the 168-bp circles in 1 mM ZnBr₂ / 10 mM MgCl₂.

Please change to caption font size

subscript "2"

The following

CCCAAAA
TTTTT

I do not like the placement of this "chart". I suggest the you make the last sentence as part of the prior paragraph and move this "chart" over to left and flush it with the copy. I also think it needs a thin rule line around it.

e-mediated ligation:

AAGGGCCAAAAGGG
CCGGTTTTTCCCGGG

High ZnBr_2 (1 mM)

DNA was imaged in either 1 mM ZnBr_2 alone or in both 1 mM ZnBr_2 and various concentrations of MgCl_2 . In all cases, regardless of the MgCl_2 concentration (which was as high as 10 mM, as shown in Figure 1), the DNA molecules were clearly kinked. The fact that zinc-induced kinking remained even when Mg^{2+} was tenfold more concentrated suggests that the binding of these two ions to DNA molecules is independent of one another. This result agrees with a well-accepted mechanism in which Zn^{2+} interacts predominantly with bases, whereas Mg^{2+} associates primarily with phosphates.

Medium ZnBr_2 (100 μM)

Figure 2 shows an AFM image of the 168-bp circles still kinked in 100 μM Zn^{2+} and 1 mM MgCl_2 , although both kinked and circular molecules have been observed under these conditions. The number of kinks was estimated to be 4. An interesting feature is the presence of “islands” around each molecule. They are likely salt clusters on which DNA molecules often deposit.

Low ZnBr_2 (50 μM)

When Zn^{2+} concentration was further reduced to 50 μM in 1 mM MgCl_2 , very few kinked molecules remained and almost all were circular (Figure 3). This result is similar to that of DNA in Mg^{2+} alone.

Based on the ~~above~~ experiments, we conclude the following.

The threshold

As there was a clear conversion from kinked to circular DNA when the Zn^{2+} concentration was reduced from 100 to 50 μM , the threshold for zinc-induced kinking evidently lies between these two values. Previous solution studies suggested that the effects of Zn^{2+} persist even to lower concentrations (10 μM) than observed here. Since Zn^{2+} must interact somewhat with the mica surface as well as DNA; however, it is possible that the activity of the Zn^{2+} is reduced substantially near the surface. It is safe to say that our threshold estimate defines an upper limit value.

Kinking mechanism

The fact that images from the same original preparation of circular DNA samples show both kinked and non-kinked structures suggests that the observed kinks cannot be attributed to nicks arising from incomplete ligation. Furthermore, when Zn^{2+} concentration varied, we did not observe any intermediate state in which the number of kinks per molecule lies between 1 and 4. Therefore, the zinc-induced kinking on a DNA molecule appears to be a cooperative process. The 4-kinks-per-molecule state is probably much more energetically favored than that of 2 or 3 kinks in this particular molecule.

Image contrast

We found an interesting correlation between the amount of salt and the apparent contrast of the DNA. As ion concentration varied from Figure 1 to Figure 3, the image contrast improved from only 0.2 nm in Figure 1 to about 1.5 nm in Figure 3. These results suggest a significant role of salt in AFM image contrast of DNA samples.

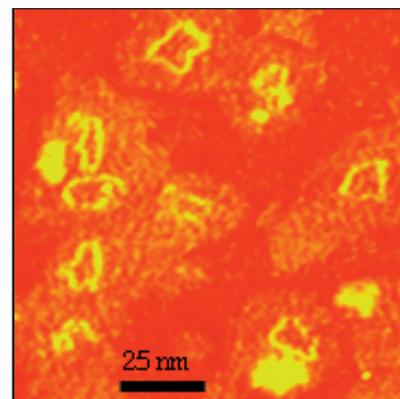


Figure 2. MAC Mode AFM image of the 168-bp circles in 100 μM ZnBr_2 / 1 mM MgCl_2 .

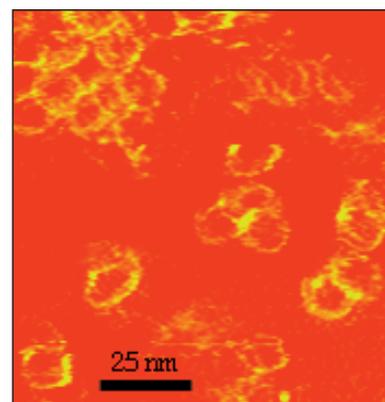


Figure 3. MAC Mode AFM image of the 168-bp circles in 50 μM ZnBr_2 / 1 mM MgCl_2 .

Summary

Using a unique AFM technology, MAC Mode, we characterized zinc-induced DNA kinking under different physiological conditions. For the first time, we were able to explore the ionic strength threshold of the kinking in a fairly direct fashion. We also showed that MAC Mode combines a major improvement in image resolution with an ability to image DNA that is relatively gently bound to the substrate. Future experiments in which a fluid-flow system is added to the MAC Mode AFM in order to obtain “online” changes of ionic environment could allow real-time kinetic studies of the kinking process.

References

1. P.K. Hansma, J.P. Cleveland, M. Radmacher, D.A. Walters, P.E. Hillner, M. Bezanilla, M. Fritz, D. Vie, H.G. Hansma, C.B. Prater, J. Massie, L. Fukunaga, J. Gurley, & V. Elings (1994) Appl. Phys. Lett. 64, 1738 -1740.
2. C.A.J. Putman, K.O.V. Werf, B.G. deGroot, N.F.V. Hulst, & J. Greve (1994) Appl. Phys. Lett. 64, 2454-2456.
3. W. Han, S.M. Lindsay, & T. Jing (1996) Appl. Phys. Lett. 69, 4111-4114.
4. W. Han, S.M. Lindsay, M. Dlakic, & R.E. Harrington (1997) Nature 386, 563.
5. J. Nickol & D.C. Rau (1992) J. Mol. Biol. 228, 1115-1123.
6. H. Sigel (1993) Chem. Soc. Rev. 22, 255-267.

This application note was created by Y. Judy Zhu, Mensur Dlakic, and Rodney E. Harrington from Arizona State University in collaboration with Keysight Technologies.

AFM Instrumentation from Keysight Technologies

Keysight Technologies offers high-precision, modular AFM solutions for research, industry, and education. Exceptional worldwide support is provided by experienced application scientists and technical service personnel. Keysight's leading-edge R&D laboratories are dedicated to the timely introduction and optimization of innovative and easy-to-use AFM technologies.

www.keysight.com/find/AFM

For more information on Keysight Technologies' products, applications or services, please contact your local Keysight office. The complete list is available at: www.keysight.com/find/contactus

Americas

Canada	(877) 894 4414
Brazil	55 11 3351 7010
Mexico	001 800 254 2440
United States	(800) 829 4444

Asia Pacific

Australia	1 800 629 485
China	800 810 0189
Hong Kong	800 938 693
India	1 800 112 929
Japan	0120 (421) 345
Korea	080 769 0800
Malaysia	1 800 888 848
Singapore	1 800 375 8100
Taiwan	0800 047 866
Other AP Countries	(65) 6375 8100

Europe & Middle East

Austria	0800 001122
Belgium	0800 58580
Finland	0800 523252
France	0805 980333
Germany	0800 6270999
Ireland	1800 832700
Israel	1 809 343051
Italy	800 599100
Luxembourg	+32 800 58580
Netherlands	0800 0233200
Russia	8800 5009286
Spain	800 000154
Sweden	0200 882255
Switzerland	0800 805353
	Opt. 1 (DE)
	Opt. 2 (FR)
	Opt. 3 (IT)
United Kingdom	0800 0260637

For other unlisted countries:

www.keysight.com/find/contactus

(BP-09-23-14)

This information is subject to change without notice.

© Keysight Technologies, 2007 – 2014

Published in USA, July 31, 2014

5989-6275EN

www.keysight.com