

# Keysight Technologies

## MAC Mode<sup>®</sup> Atomic Force Microscope for Precision Interfacial Force Measurements

Application Note

## Introduction

Since its introduction, it has been recognized that the atomic force microscope (AFM) offers a unique ability to probe surface forces with nanometer resolution and that these surface forces play a vital role in the operation of the AFM (1). This application note describes a new approach to AFM in which a force is applied directly to an AFM tip using a magnetic field (MAC Mode AFM), rather than by pressing the tip into the sample with the scanning transducer. This method eliminates the spurious resonances seen in fluid tapping mode (3), leading to substantial improvements in image quality (6). Here, we discuss the application of MAC Mode AFM to the measurement of surface forces.

Ordering of the liquid at the interface with a solid is a phenomenon of fundamental importance in fields such as tribology and electrochemistry and it has been studied extensively with a surface force apparatus (SFA) (7). This liquid structuring can also have a profound effect upon the operation of an AFM, causing the tip to hop between different stable configurations in the ordered region of the fluid (8). The AFM has been used to study structure in interfacial water (9) and organic liquids (5, 10). In contrast to the SFA, the AFM suffers from a lack of reproducibility of one surface (the tip), but it has a number of advantages, among them the ability to produce high-resolution images of the surface along with high vertical (z) sensitivity and low mechanical drift. This application note illustrates how highly reproducible and accurate data can be obtained despite variability in AFM tips.

If the microscope is operated with a probe that is oscillated at an amplitude  $x$  far from the surface, then, provided the oscillation frequency is far below mechanical resonances, the surface stiffness at some distance  $z$  from the surface,  $S(z)$ , is given by  $S(z) = k [ x_0 / x(z) - 1 ]$  where  $k$  is the spring constant of the cantilever and  $x(z)$  is the amplitude of the tip oscillation (5). A number of simple models permit extraction of quantities such as Young's modulus from the surface stiffness (1, 11).

A substantial signal-to-noise advantage is obtained for A/C mode AFM in fluids if the tip is oscillated by direct application of a magnetic force, as opposed to indirect mechanical excitation with an acoustic transducer. O'Shea et al. (5) used a magnetically-oscillated-tip microscope to show evidence of up to seven ordered layers of liquid at an interface with graphite. However, the spacing between layers of the quasi-spherical molecule studied, octamethylcyclotetrasiloxane (OMCTS), was substantially smaller (5 Å) than the known smallest dimension of the molecule (8 Å); thus, an orientation could not be determined. We have repeated these measurements with a microscope designed for stable operation in fluids. For both OMCTS and another small molecule (mesitylene) we obtained data in excellent agreement **Subscript highlighted text**

We chose OMCTS ( $(\text{CH}_3)_2\text{SiO}_4$ ) and mesitylene ( $(\text{CH}_3)_3\text{C}_6\text{H}_3$ ) because they are relatively rigid and their dimensions are well characterized. The oblate spherical model of OMCTS has a major diameter of 1.0 to 1.1 nm and a minor diameter of 0.7 to 0.8 nm (12). Ordered layers have been detected on mica with a period of 0.8 to 1.0 nm using the SFA (13). Mesitylene or 1,3,5-trimethylbenzene molecule has a benzene ring with each -H at the 1<sup>st</sup>, 3<sup>rd</sup>, and 5<sup>th</sup> carbon locations replaced by a -CH<sub>3</sub> group. Crystal structure data show that the ring-to-ring spacing in stacked layers is 3.6 Å in the solid, and that the end-to-end length of the molecule along the major axis is a little over 6 Å. OMCTS was obtained from Sigma (St. Louis, Missouri); mesitylene was obtained from Alfa Aesar (Ward Hill, Massachusetts). Highly oriented pyrolytic graphite (HOPG) was obtained from Advanced Ceramics (Cleveland, Ohio).

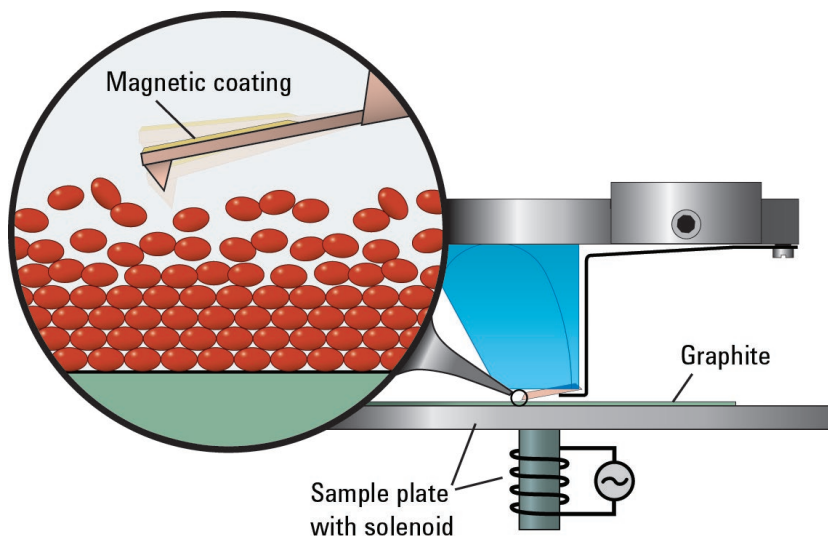



Figure 1. Schematic setup of the MAC Mode AFM.

We used a Keysight AFM operated in MAC Mode (see Figure 1). The spring constant of the MAClevers (magnetically coated cantilevers) used was typically 0.5 nN/nm. Atomic force spectra (AFS) were obtained by driving the force input of the controller with the demodulated amplitude signal from the MAC Mode and setting the controller to record force curves. The tip oscillation frequency was between 200 to 500 Hz, far below the 28 kHz resonance of the cantilever in fluid.


The fluid cell and cantilever holder were thoroughly cleaned before each experiment. OMCTS was always stored in a desiccator filled with  $N_2$  at 4 °C. It was warmed to ambient temperature before being taken out. In each experiment, graphite was cleaved immediately before being covered by fluid. All experiments were carried out at room temperature in air. The characteristics of AFS data presented here did not change over hours of measurements in our experiments.

Figures 2(a) and (b) are typical force spectra from two independent experiments at OMCTS-graphite interfaces. They demonstrate the reproducibility of the spectra obtained with different tips. Both approaching and retracting curves were recorded. In Figure 2(a), seven layers of OMCTS molecules at the interface are clearly detected. The average spacing between two adjacent layers, which reflects the size of an OMCTS molecule along the direction perpendicular to the layers, was measured to be 8.7 Å for approaching and 7.9 Å for retracting. The difference was due to thermal drift along the approaching direction. In Figure 2(b), seven molecular layers are also detected. Thermal drift was less in this case. The average spacing between two adjacent layers was 8.2 Å for approaching and 8.0 Å for retracting. After analyzing 22 force spectra from three separate experiments, we found that the average size of OMCTS molecules along the direction of the ordered layers at the interface was  $8.2 \pm 0.3$  Å, consistent with the minor axis of the molecule according to the oblate spherical model of OMCTS (12) and considerably smaller than the major axis. These results show that the layers stack along the minor axis of the molecules.


Figures 3(a) and (b) show two measurements of typical force spectra at mesitylene-graphite interface. Once again, the general features of the AFS are seen to be independent of the tip used. Figure 3(a) shows eight molecular layers on approach and six on retraction, presumably because of disruption of the fluid structure by the tip. The average period of the layers was measured to be 4.2 Å for approaching and 4.3 Å for retracting. Figure 3(b) shows less layers (six approaching and five retracting) and slightly

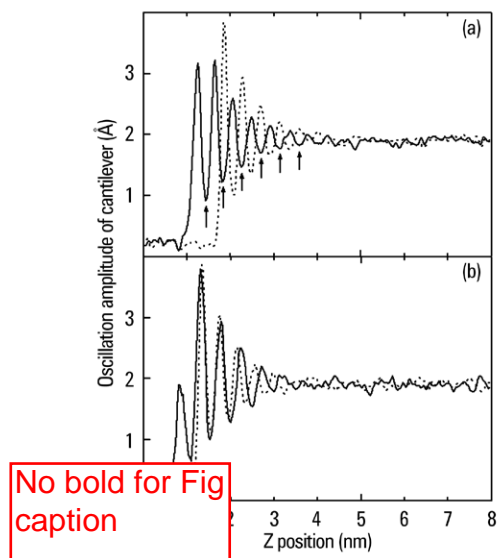

**Figure 2.** Two independent measurements of force spectra of the MAC Mode AFM at OMCTS-graphite interface. The amplitude of oscillation of the cantilever oscillates in both approaching (solid line) and retracting (dotted line) curves in the region of a few nanometers away from the surface due to ordered layers of OMCTS molecules. The period, 8.2 Å, precisely reflects the layer thickness of OMCTS molecules. (a) Driving frequency = 500 Hz; scan rate = 2.8 nm/sec. (b) Driving frequency = 200 Hz; scan rate = 1.6 nm/sec. The arrows on this and subsequent plots correspond to the repulsive-force maxima.

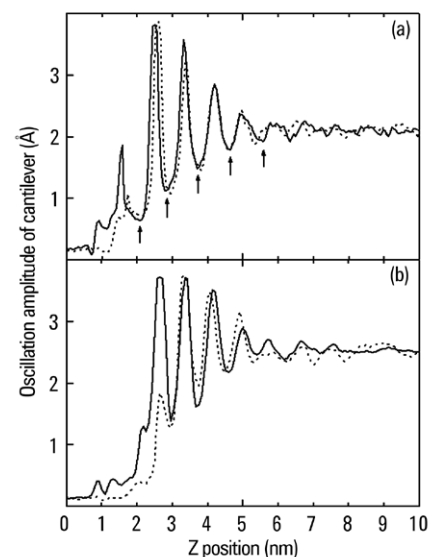

**Figure 3.** Two typical force spectra at mesitylene-graphite interface in both approaching (solid line) and retracting (dotted line) processes. Up to eight layers of mesitylene were observed. The oscillatory period, 4.4 Å, indicated that the ordered mesitylene rings arranged parallel to the graphite surface. For both (a) and (b): driving frequency = 500 Hz; scan rate = 2.4 nm/sec.

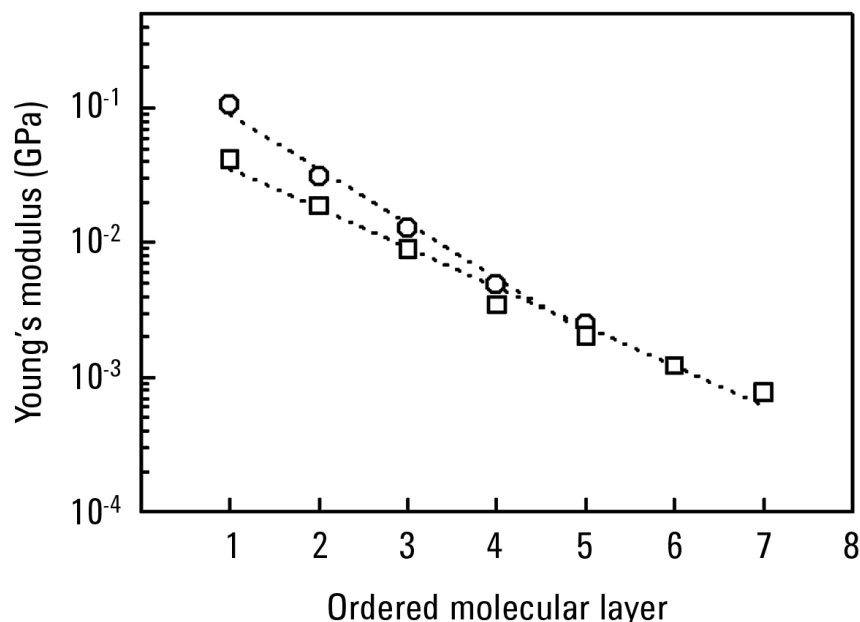

**Figure 4.** The Young's modulus of OMCTS (circles) and mesitylene (squares) at centers of ordered layers. The lines are best-fit lines, indicating exponential decay of the Young's modulus with the number of molecular layers. Data from Figure 2(a) and Figure 3(a).

higher thermal drift. The average periods for approaching and retracting were 4.7 and 4.1 Å, respectively. After averaging 22 force spectra in three separate experiments, we found the size of mesitylene molecules along the direction perpendicular to the ordered layers is  $4.5 \pm 0.2$  Å. Given the crystal structure data (14), these results suggest that the mesitylene rings lie flat with respect to graphite surface.

The liquid density profile at liquid-solid interface oscillates due to the constraint of the solid surface and interactions between the surface and the liquid (15), thus resulting in a corresponding oscillation in the surface stiffness. The amplitude of oscillation of the tip was set at about 2 Å far from the surface. It can be seen that the layering introduces regions of both positive stiffness (decreased amplitude) and negative stiffness (increased amplitude).

For discrete molecular layers, it is inappropriate to map the liquid Young's modulus (compliance) continuously away from the surface because negative stiffness occurs around centers between two neighboring layers. However, it is possible to estimate the Young's modulus at centers of molecular layers where the tip sensed a local maximum repelling force. These points are marked with arrows on both Figures 2 and 3. Note that the amplitude of these minima changes in a consistent way from run to run, while the peak positions, corresponding to negative stiffness, do not.

The minima correspond to stable repulsive mechanical contact <sup>superscript "th"</sup> with simple model-based Hertz elasticity (11). An effective Young's modulus for the  $n$ th layer,  $E_n^*$ , can be obtained from the surface stiffness for the  $n$ th layer,  $S_n$ , according to  $E_n^* = S_n [S_n / 6RL_n]^{1/2}$  where  $R$  is the radius of the tip and  $L_n$  is the load on the  $n$ th layer (11).  $L_n$  is calculated from  $L_n = x_n k$  where  $x_n$  is the amplitude decrease for the  $n$ th layer. Assuming that the radius of the tip is 30 nm and its spring constant is 0.5 nN/nm, the Young's modulus at the center of each ordered molecular layer was plotted in Figure 4 for both OMCTS (circles) and mesitylene (squares). The best-fit lines yield decay lengths of 1.07 layers for OMCTS, corresponding to 8.8 Å, and 1.48 layers for mesitylene, corresponding to 6.5 Å.

Please move Keysight to the next line to be with Technologies

In conclusion, we have used AFS of the MAC Mode AFM to probe ordered liquid molecules of OMCTS and mesitylene at the graphite surface with sub-angstrom vertical resolution. Elastic compliance at the center of each layer was found to exponentially decrease versus the number of layers. Combined with the high lateral resolution feature of the MAC Mode AFM (6), we believe this study will open a possible way of precisely identifying properties of molecules at surfaces and interfaces from solid to liquid.

## References

1. J.B. Pethica and W.C. Oliver, *Physica Scripta* **T19**, 61 (1987).
2. S.M. Lindsay, Y.L. Lyubchenko, N.J. Tao, Y.Q. Li, P.I. Oden, J.A. DeRose, and J. Pan, *J. Vac. Sci. Technol.* **11**, 808-815 (1993).
3. Wenhai Han, S.M. Lindsay, and Tianwei Jing, *Appl. Phys. Lett.* **69**, 4111 (1996).
4. M.A. Lantz, S.J. O'Shea, and M.E. Welland, *Appl. Phys. Lett.* **65**, 409-411 (1994).
5. S.J. O'Shea, M.E. Welland, and J.B. Pethica, *Chem. Phys. Lett.* **223**, 336 (1994).
6. Wenhai Han, S.M. Lindsay, M. Dlakic, and R.E. Harrington, *Nature* **386**, 563 (1997).
7. J.N. Israelachvili, *Intermolecular and surface forces*, Academic Press, New York, 1992.
8. K. Koga and X.C. Zeng, *Phys. Rev. Letts.* **79**, 853 (1997).
9. J.H. Hoh, J.P. Cleveland, C.B. Prater, J.-P. Revel, and P.K. Hansma, *J. Am. Chem. Soc.* **114**, 4917 (1992).
10. S.J. O'Shea, M.E. Welland, and T. Rayment, *Appl. Phys. Lett.* **60**, 2356 (1992).
11. R.M. Overney, E. Meyer, J. Frommer, H.-J. Güntherodt, M. Fujihira, H. Takano, and Y. Gotoh, *Langmuir* **10**, 1281 (1994).
12. D.W. Scott, *J. Am. Chem. Soc.* **68**, 2294 (1946).
13. R.G. Horn and J.N. Israelachvili, *J. Chem. Phys.* **75**, 1400 (1981).
14. Tor Dahl, *Acta. Chem. Scand.* **25**, 1031 (1971).
15. F.F. Abraham, *J. Chem. Phys.* **68**, 3713 (1978).

### AFM Instrumentation from Keysight Technologies

Keysight Technologies offers high-precision, modular AFM solutions for research, industry, and education. Exceptional worldwide support is provided by experienced application scientists and technical service personnel. Keysight's leading-edge R&D laboratories are dedicated to the timely introduction and optimization of innovative and easy-to-use AFM technologies.

[www.keysight.com/find/AFM](http://www.keysight.com/find/AFM)

For more information on Keysight Technologies' products, applications or services, please contact your local Keysight office. The complete list is available at: [www.keysight.com/find/contactus](http://www.keysight.com/find/contactus)

#### Americas

Canada	(877) 894 4414
Brazil	55 11 3351 7010
Mexico	001 800 254 2440
United States	(800) 829 4444

#### Asia Pacific

Australia	1 800 629 485
China	800 810 0189
Hong Kong	800 938 693
India	1 800 112 929
Japan	0120 (421) 345
Korea	080 769 0800
Malaysia	1 800 888 848
Singapore	1 800 375 8100
Taiwan	0800 047 866
Other AP Countries	(65) 6375 8100

#### Europe & Middle East

Austria	0800 001122
Belgium	0800 58580
Finland	0800 523252
France	0805 980333
Germany	0800 6270999
Ireland	1800 832700
Israel	1 809 343051
Italy	800 599100
Luxembourg	+32 800 58580
Netherlands	0800 0233200
Russia	8800 5009286
Spain	800 000154
Sweden	0200 882255
Switzerland	0800 805353
	Opt. 1 (DE)
	Opt. 2 (FR)
	Opt. 3 (IT)
United Kingdom	0800 0260637

For other unlisted countries:

[www.keysight.com/find/contactus](http://www.keysight.com/find/contactus)

(BP-09-23-14)

This information is subject to change without notice.

© Keysight Technologies, 2007 - 2014

Published in USA, July 31, 2014

5989-6377EN

[www.keysight.com](http://www.keysight.com)